




# VIAL2BAG ADVANCED®

UNIVERSAL COMPATIBILITY WITH ALL MANUFACTURERS' IV BAGS IN 50ML, 100ML, OR 250ML THAT CONTAIN AN ISO STANDARD IV PORT\*

The Vial2Bag Advanced® 20mm admixture device enables reconstitution and transfer of a drug between a vial and an IV bag prior to administration to the patient. The innovative needle-free and easy-to-use design is ideal for immediate use and can optimize pharmacy clean room time, while assisting in standardizing admixture procedures to maximize efficiency of point-of-care nursing.

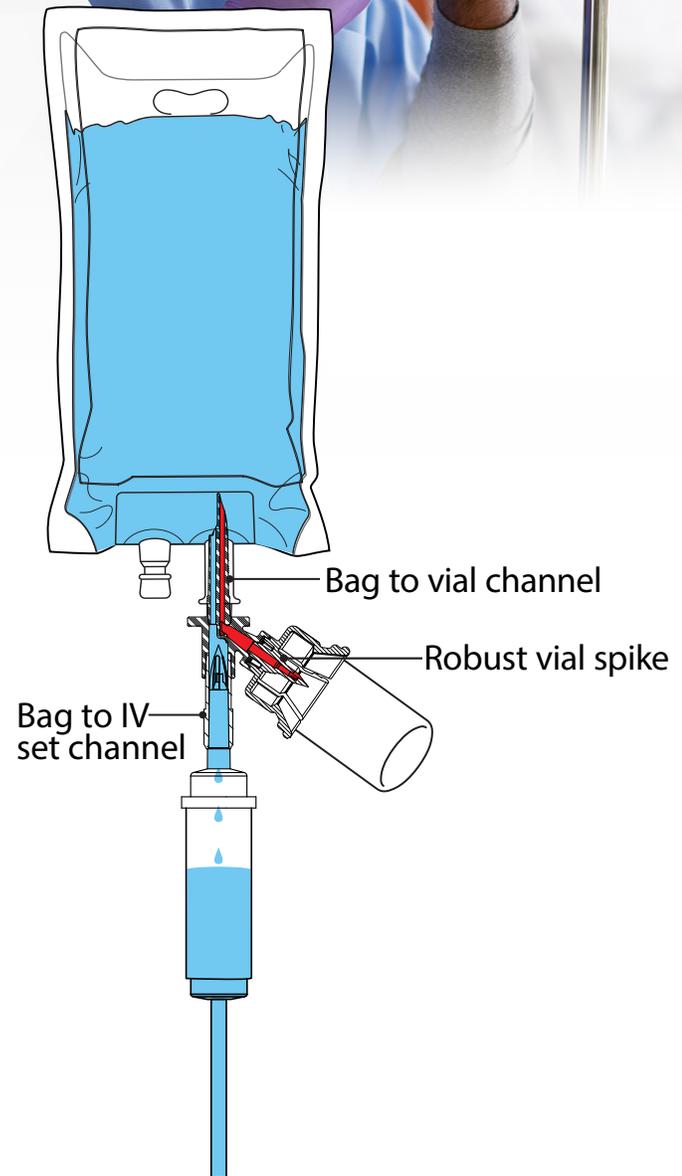

## NEW DESIGN

- Dual channel design provides dedicated fluid pathways into and out of the IV bag
- Enhanced vial spike design improves the connection between vial and Vial2Bag Advanced® 20mm admixture device adapter and the transfer of medication†
- Needle-free device helps to reduce potential needlestick injuries

\*ISO 8536-4 standard IV spike  
†As Compared to Vial2Bag® DC



The Vial2Bag Advanced® admixture device can be an effective tool to assist healthcare providers in immediate use preparation, in accordance with USP <797>.

### Optimized Patient Care

- Provides a simple solution for just-in-time IV medication preparation, improving timeliness at point of care
- Easy-to-implement immediate use solution

### Effective Waste Reduction

- Can reduce medication and IV fluid waste
- Universal compatibility with all manufacturers' IV bags in 50mL, 100mL, or 250mL that contain an ISO standard IV port
- Suitable for use with liquid and powder medications with standard 20mm closure vials<sup>‡</sup>

<sup>‡</sup>Note: The Vial2Bag Advanced® 20mm admixture device is 510(k) cleared by the United States Food and Drug Administration (FDA). The use of the Vial2Bag Advanced® 20mm admixture device should not be interpreted as modifying, extending, or superseding a drug manufacturer's labeling recommendations for storage and expiration dating, unless otherwise limited by USP <797> compounding standards. Refer to drug manufacturer's labeling and use instructions for recommendations, USP <797>, and applicable institution policy for shelf life and sterility information of reconstituted product and admixture device compatibility. Compatibility of the Vial2Bag Advanced® 20mm admixture device with all drug products has not been confirmed. Do not use the Vial2Bag Advanced® 20mm admixture device with lipids. Failure to follow the instructions provided may result in inadequate medication reconstitution, dilution, and/or transfer, possibly leading to overdose or underdose and/or delay in therapy. Products shown are for INFORMATION purposes only and may not be approved for marketing in specific regions. Please contact your West Pharmaceutical Services, Inc. (West) representative for product availability.

**Rx Only**

## Distribution and Ordering



**Vial2Bag Advanced®  
20mm Admixture Device**  
Item # 36098150

**Quantity**  
50/Box  
150/Case



**Exclusively Distributed by: Progressive Medical, Inc.**  
997 Horan Dr. • Fenton, MO 63026  
**Ordering Information:**  
Phone: 314-961-5786 • Fax: 314-961-4535  
cs@progressivemedinc.com



## INNOVATION IN RECONSTITUTION AND TRANSFER SYSTEMS

Scan to learn more about our products and services

For more information visit [www.westpharma.com](http://www.westpharma.com)